

B L U M I L - G O

BLACK
LINES!



Blumil -GO is extremely easy to use. You stay seated in your own wheelchair, approach the Blumil -GO, press one button and you are connected and ready to go. Blumil -GO will convert your manual wheelchair into an electric handbike in a matter of seconds. **It's that simple.**



DESIGN

Design across our full product range is of paramount importance to us. The Blumil -GO is testament to this! We wanted to design a product that was not only practical, but one that looks slick, compact and trendy. No unnecessary cables, attachments and parts.






P **PATENTED**


PATENTED CONNECTING MECHANISM

The Blumil -GO has one incredibly unique feature that has never been seen before - patented connecting mechanism. We decided to do things a little different in order to make a HUGE difference to the user. There is no need to attach any additional parts to your own wheelchair. No troublesome procedure of having to connect to the frame of your wheelchair by lifting the whole construction. With the Blumil -GO you just press a button and an electric motor does the rest for you. 5 seconds later you are ready to GO with absolutely no assistance from anyone. You can be fully independent, no matter how much strength you have in your hands.

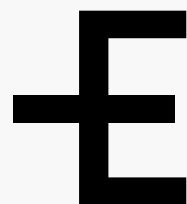


EASY TO MANEUVER

With a turning circle of only 1,8m, the Blumil -GO is extremely easy to maneuver. This is a particularly important feature indoors as well as outdoors. It is amazingly easy to turn in tight spaces and small rooms. Turning around when you are driving on a path requires minimal space.

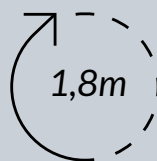

Turning Circle



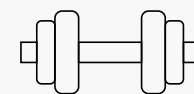

13kg - *Weight*



Travel Easy

The possibility of easy transport is a key feature for Blumil -GO users. Because of this, we ensured the Blumil -GO would have a really compact size, with the ability to fold effortlessly with the touch of a button.


Take it anywhere with you.

LIGHT AND COMPACT

With a total weight of 13 kg, INCLUDING the battery, the Blumil -GO is practical enough to travel with you, whether by car or plane. There are NO additional parts that you need to install to your wheelchair. The complete system that connects to your wheelchair in 5 seconds is that light.

TRAVEL BY PLANE

It is especially great for traveling. It folds easily, fits in a custom dedicated bag and you can take it with you wherever you go!

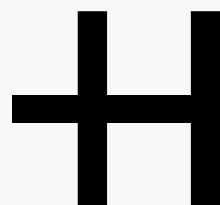


HANDLING OBSTACLES

BLUMIL-GO



With the Blumil -GO you move on 3 wheels (2 rear wheels of your wheelchair and one front wheel of the Blumil -GO). This makes the ride over any surface much more comfortable compared to moving on 4 wheelchair wheels, as the casters (small front wheels) are slightly raised.

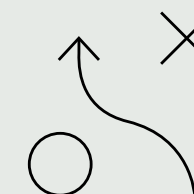



D

DOUBLE BRAKING


Blumil -GO has a unique double braking system.

An easy to operate electromagnetic brake for people with limited hand movement and a second powerful main brake operated by a brake lever.



GREAT FOR
QUADRIPLLEGICS

Blumil -GO is easy to use, drive, connect and set up. For 99% of users, the Blumil -GO can be unboxed, and set up within 15 minutes with no adaptations to the wheelchair. Starting from the unique patented connecting system, which is activated by just a touch of a button to driving.

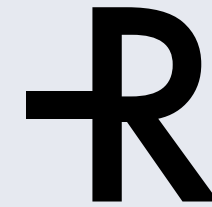


REMOVABLE BATTERY

Thanks to its reliable 280 Wh LG battery pack, the Blumil -GO has an impressive range of approximately 30 km per charge. You can easily remove it and take it with you. You can charge it on the go and it additionally has a USB outlet so you can charge all your mobile devices.

You have the biggest power bank ever!



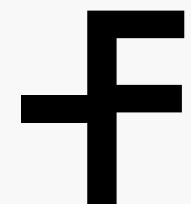

FEATURES



SPEED

Blumil -GO has 3 speed settings. Choose whichever you feel comfortable with. The standard setting has a top speed of 6 km/h, which is great when you are starting to learn how to drive or if you want to travel at walking pace with friends and family. The middle setting increases the maximum speed to 16 km/h. This allows the user to swiftly commute around city centers and parks alike and is ideal for day to day use. For those that love an adrenaline rush there is a TURBO setting which increases the Blumil -GO's speed to a maximum of 30 km/h allowing you to keep up the pace with any bike out there, just be careful.

This is for expert drivers!



FITS RIGID FRAME WHEELCHAIRS


Blumil -GO was specifically designed for rigid frame wheelchairs. It should easily fit all brands and models of rigid frame wheelchairs from the most popular manufacturers like: Panthera, Kuschal, Invacare, Meyra, Progeo and more.

Please note that it will not fit folding wheelchairs.



TAKE IT OUT OF THE BOX AND GO

The Blumil -GO is not only easy to use, but also amazingly easy to set up to fit your wheelchair.

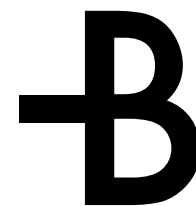


BUILT IN LIGHTS

For added safety, the Blumil -GO has two sets of lights.

The ambient light turns on automatically as you turn the Blumil -GO on, instantly making you visible to others.

Additionally, you can switch on the main headlights, when you need extra light at night to guide your path.

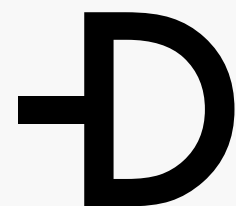

DIMENSIONS



47cm

92cm




88cm



50cm

110cm

Three speed modes	6 km/h, 16 km/h, 30 km/h
Range	25 km
Removable battery	Li-ion 280 Wh
Total weight	13 kg
Wheel diameter	8.5"
Maximum load	Max. 100 kg
Lights	LED
Tire	Pneumatic tubeless puncture proof
Connecting mechanism	fully electric
Certification	CE, medical device in EU
Climbing capability	Up to 10 %
Digital LCD-display	current speed and speed mode, charge level, total distance
Dual braking system	drum brake and electromagnetic brake

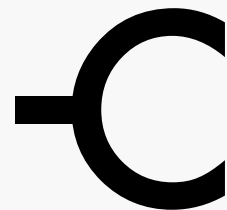

CERTIFICATION

B L U M I L -GO



CERTIFICATION

Blumil -GO is a certified product. As a medical device it meets all the requirements, regulated under EU law. The CE mark itself is a determinant of certain safety and quality standards and provides information for all consumers and proof that we meet all strict requirements.



MADE IN THE EU

Blumil -GO is of very high quality. It is manufactured in the EU to the highest standards.



B L U M I L -GO

CONTACT US

Office address:

ul. Geldowa 4B/11, 01-211 Warsaw, Poland

Email:

info[@]blumil.com

www.blumil.com